



# Safe Electrosurgery


How to increase your safety and thus of the patient?

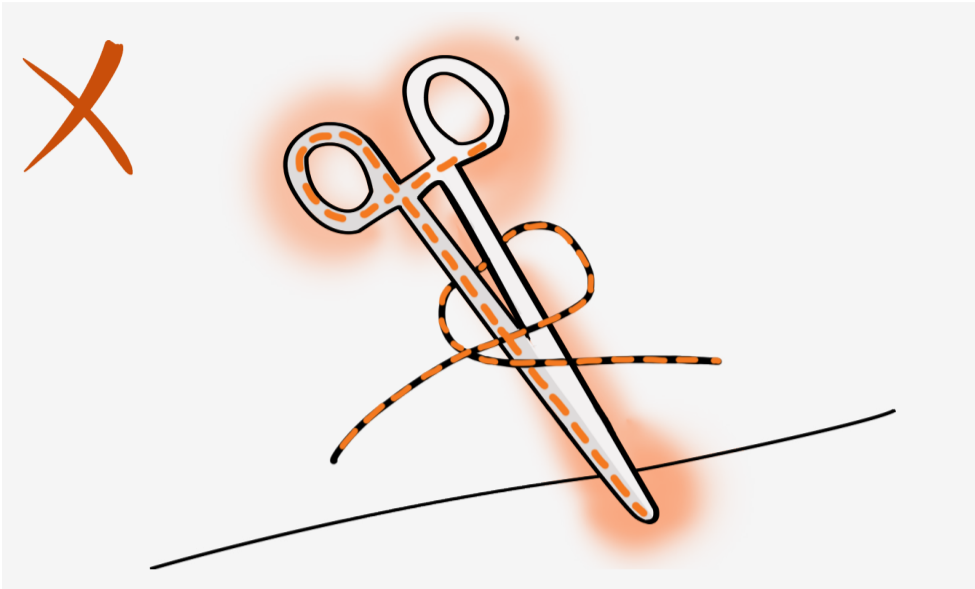

**1** Do not wrap the handheld cord around any metal instruments!  
This could activate these instrument and lead to patient burns.

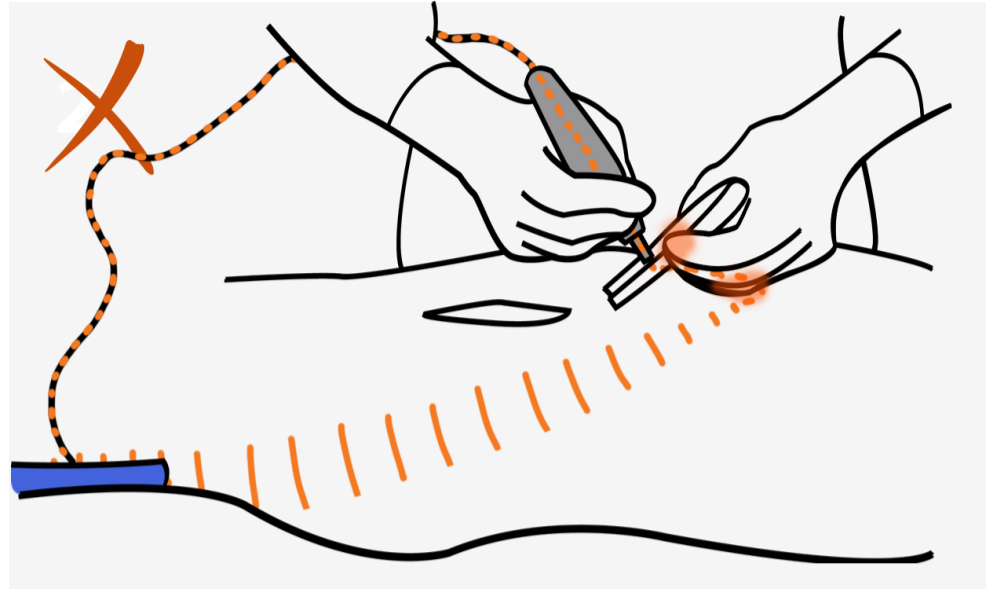


**2** Do not activate the system when the active electrode is not in contact with tissue!  
This can lead to glove burns (operator and patient burns).

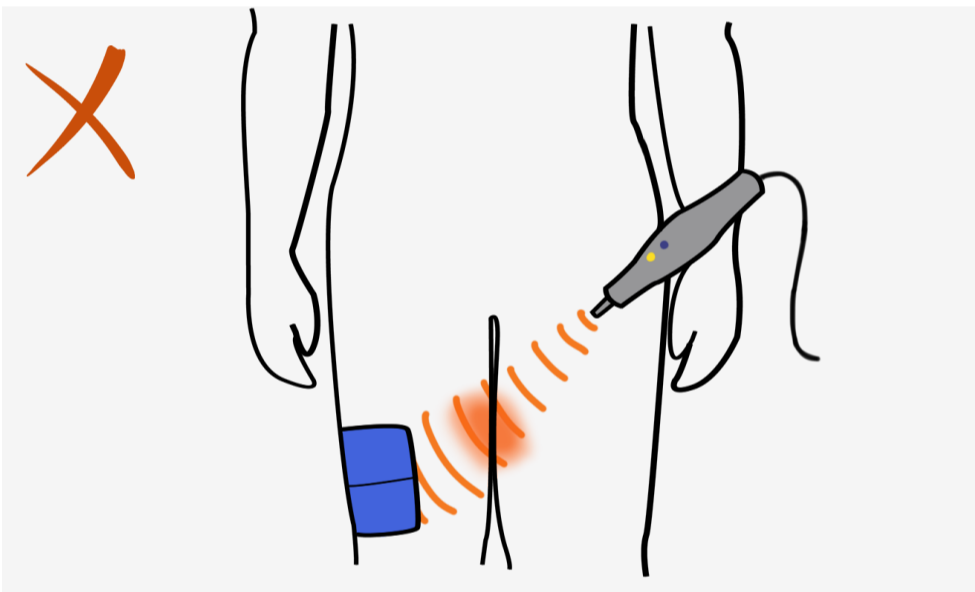


**3** Avoid skin-to-skin contact.  
This could lead to patient burns.



**4** Make sure that the patient is not in contact with grounded metal objects!  
This could lead to patient burns.

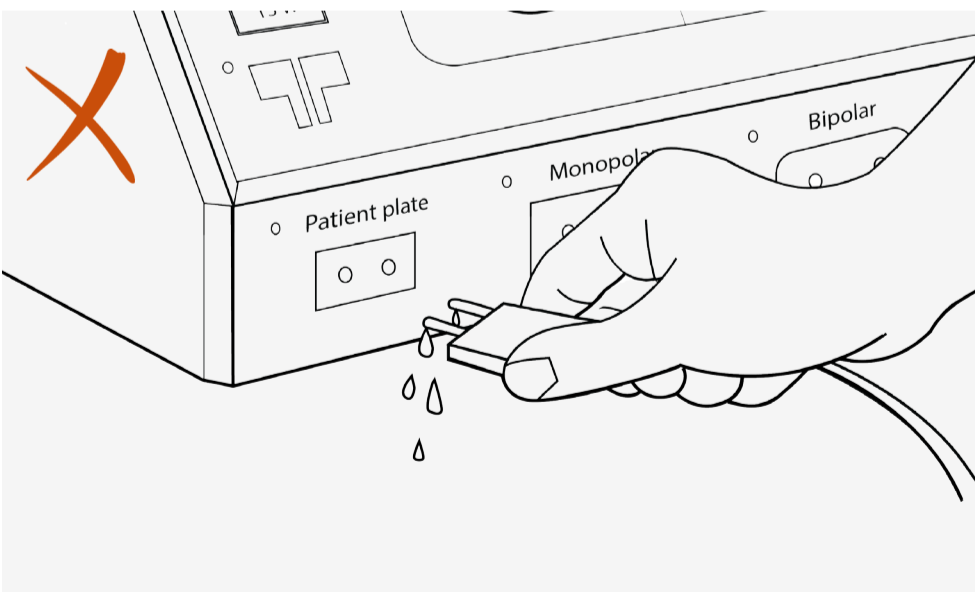

**5** Do not plug wet accessories into the outlets of the HF generator.  
This could lead to electrocution!

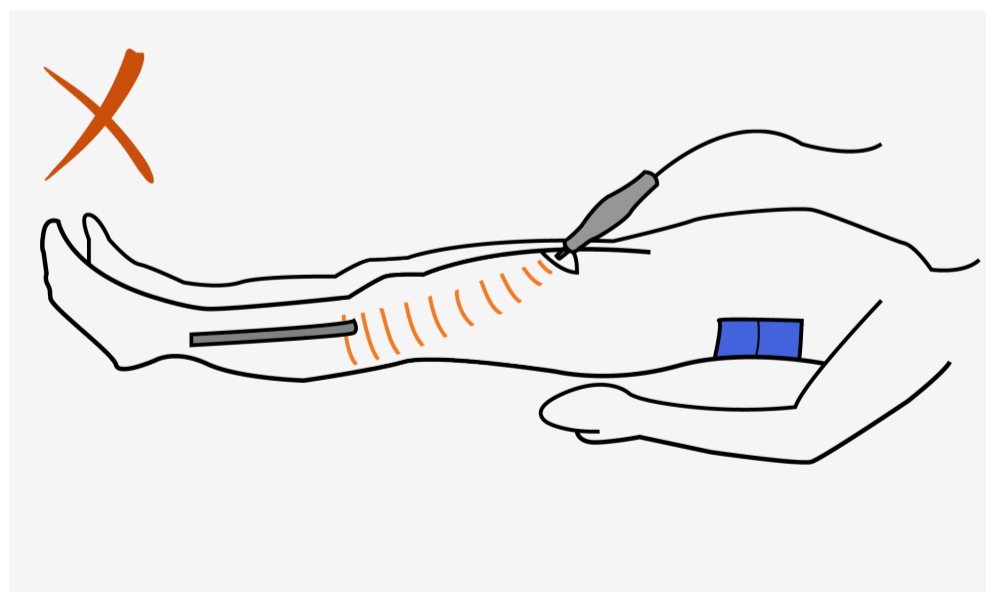

**6** Make sure that the patient has no metal objects inside his body. This could influence the current flow and lead to burn wounds.